




PRECISION LOUPES, FITTED TO YOU

See every detail Keep your posture


Deflection prism optics let you work looking straight ahead.

Look straight ahead, not down

Conventional loupes converge your eyes toward a single close point, and to see the field you tip your head down and hold it there. PENTAX Loupes are built around a different idea: your eyes stay level, and the optics do the bending.

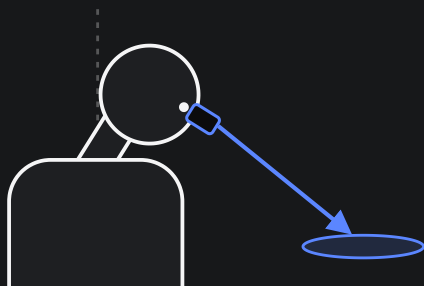

Conventional loupes. The head tips forward so the eyes can converge on one close point, and stays there for the whole procedure.

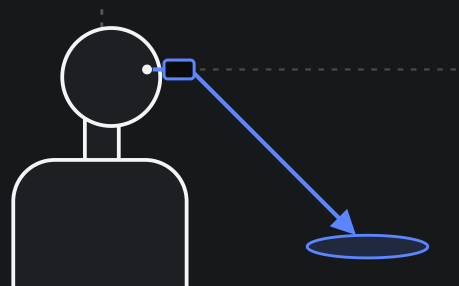

PENTAX Loupes. The head stays level and the eyes look straight ahead. Prisms inside the barrels bend the optical path down to the field.

Deflection prisms

Prisms inside each barrel redirect the optical path, much like working with a microscope. You look straight ahead at your work instead of angling your eyes down, which helps you hold a natural, upright posture through long procedures.

Parallel viewing

Conventional loupes converge your eyes toward a single close point. PENTAX Loupes are built for parallel viewing instead, designed to reduce eye convergence and help support an upright working posture.



"With PENTAX loupes, the posture of my neck really changed. I can look straight ahead, just like working with a microscope, and I enjoy parallel viewing with a very large field of view."

Dr Robotti · President, The Rhinoplasty Society of Europe

Built around you

Every measurement that matters is adjustable, and it is set for you at your fitting.



25–70 cm

Working distance

Set to match your height, your chair and the way you sit or stand at work.



52–78 mm

Pupillary distance

Matched to your own eyes, not a factory average.



Adjustable line of sight

A neutral posture

Working distance and where your eyes sit to look through the lenses are completely adjustable, so your head, neck and back stay in a neutral position.



Optional lenses

Protective lenses

Protective lenses fit inside the frame. Choose them at your fitting, or use the frame on its own.

Adjustable later, too. Pupillary distance and working distance can be reset as your eyesight or working posture changes, without returning the loupes to the factory.

Comfortable through long procedures

The frame, the optics and the flip-up barrels are designed to stay comfortable through long procedures.




Balance point

A lightweight frame with balanced weight distribution, so the loupes sit evenly and stay comfortable through long procedures rather than dragging on your nose.



Conventional
loupes

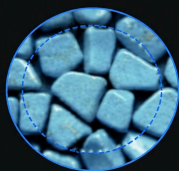

PRISMVUE

A wide field of view

Japanese optics give a sharp, distortion-free image across a wide field, with the depth perception fine work needs. At 2.5× the field of view is 141 mm.



Flip up, flip down

Flip the barrels out of the way to switch instantly between magnified and normal vision, without taking the frame off.



"With the PENTAX Loupes I have a wide field of vision and they are light to wear. You can keep your focus while looking straight ahead, which minimises strain on my neck."

Dr Dhooqhe · Surgeon

MAGNIFICATION

Choose your magnification

Five fixed magnifications, each with its own field of view. Your specialist helps you choose at your fitting.

MAGNIFICATION	FIELD OF VIEW	GUIDANCE
2.5×	141 mm	For those who favour field of view and depth of vision: general surgeons, urologists and hygienists.
3.0×	120 mm	For general dentists and ENT surgeons.
3.5×	106 mm	For anyone who wants a good middle ground between field of view and magnification.
4.0×	98 mm	For those who need good attention to detail while maintaining a decent field of view.
5.0×	80 mm	For those working on very small, fine details, such as endodontics or microsurgery.



*These are estimates, results per person may vary.

2.5×–3.5×

Dentistry

Restorative and hygiene work at 2.5× to 3.5×. Aesthetic dentists often go to 5.0×.

3.5×–5.0×

Surgery

Fine, delicate procedures at 3.5× and above, across surgical specialities.

2.5×–3.5×

Veterinary

General and referral practice, depending on procedure and patient size.

Made for your eyes

The frame adjusts to your vision. Wear the loupes on their own, or over the glasses you already have.



Worn directly

Many clinicians do not need separate glasses at all. The frame adjusts directly to your vision.



Over your own glasses

A different nose pad makes room for your existing glasses underneath.

Light where you are looking

An optional LED headlight clips to the loupes: **wired** for constant power from a beltpack, or **wireless** for complete freedom of movement. Brackets are available for other LED brands.

COLOUR

Four colours, pick yours



Silvergold



Blue



Black



Red

UK SERVICE

From enquiry to fitting

01

Tell us what you need

Reply to this email, or build your loupes online. It takes about a minute.

02

Talk it through

A specialist contacts you to discuss magnification, working distance and fit.

03

Demonstration & fitting

We arrange a demonstration and fitting with you in the UK.

How PENTAX Loupes compare

Most loupes you will have seen are through-the-lens: the barrels are fixed into the lens to one set of measurements. PENTAX Loupes take a different approach.

	PENTAX Loupes	Conventional through-the-lens loupes
Pupillary distance	You adjust it, 52–78 mm	Fixed at the factory
Working distance	You adjust it, 25–70 cm	Fixed at the factory
If your eyesight changes	Re-adjust at any time, without returning them to the factory	Built to one set of measurements; changes mean going back to the maker
Line of sight	Straight ahead. Prisms bend the optical path down to the field	Eyes converge and look down through angled barrels
Normal vision mid-procedure	Flip the barrels up, out of your sightline. The frame stays on	Look over the top of the fixed barrels
Your own glasses	Wear them underneath, with a different nose pad	Prescription is ground into the carrier lens
Fitting	Demonstrated and fitted on you, in the UK	Measured, then made to order
Warranty	Lifetime warranty	Varies by maker

The through-the-lens column describes the category in general; individual makers vary. Lifetime warranty terms are provided with your quotation.

Fully adjustable

Every measurement that matters is yours to set, and to reset.

Lifetime warranty

On every pair, for as long as you wear them.

Fitted in the UK

Demonstrated on you before you decide, then set up on your face.

WATCH

See them in use

Films from PENTAX Loupes. Click a card to watch, or point your phone camera at the code.



The PENTAX Loupes film

youtu.be/dZ84SZZZHM



Setting up your PENTAX Loupes, step by step

youtu.be/2v0Diu6g2wc



A day at PENTAX Loupes headquarters

youtu.be/aAK3cu8YtAw



Dr Cameron McIntosh on his PENTAX Loupes

youtu.be/KJ066v9qihg



Dr Enrico Robotti in theatre with PENTAX Loupes

youtu.be/WZKcJVDRXQc



What clinicians ask us

Is the magnification adjustable?

No. Magnification is fixed when the loupes are built, which gives a consistent, predictable image. Pupillary distance and working distance are adjustable by you.

Can I adjust the fit later?

Yes. Pupillary distance and working distance can be reset as your eyesight or working posture changes, without returning the loupes to the factory.

Can I wear them with my own glasses?

Yes. Many clinicians do not need separate glasses at all. If you do, a different nose pad makes room for your own glasses underneath.

Which magnification suits my work?

As a guide, dentists commonly use 2.5× to 3.5× (aesthetic dentists often 5.0×), surgeons 3.5× to 5.0×, veterinary surgeons 2.5× to 3.5×. Your specialist helps you choose at your fitting.

Can I try them before I decide?

Yes. We arrange a demonstration and fitting with you in the UK, so you see the view for yourself before anything is ordered.

What do they cost?

It depends on the magnification and on options such as LED lighting. Your specialist prepares a personalised quotation for your configuration.

What warranty do they come with?

A lifetime warranty. The terms come with your quotation.

Are LED lights available?

Yes, wired or wireless. If you already own a light you like, brackets are available for other LED brands.



"These loupes are just great. I felt a difference in my posture the very first time I tried them!"

Dr Hsu · Dentist



The best way to understand PENTAX Loupes is to wear them

Request a demonstration and fitting in the UK. A specialist will talk through magnification, working distance and fit with you, and prepare a personalised quotation for your configuration.

WEBSITE

PentaxLoupes.co.uk

EMAIL

Info@PentaxLoupes.co.uk

AFTER YOUR FITTING

Every pair comes with a passport

A private page and a QR card in the case, holding your setup, your care guidance and a one-tap line to the UK team. Magnification, working distance and pupillary distance are recorded as they are set at your fitting.

Scan the card with any phone. No app, no account, no password. Use it for a quick check-in, a support request, or to invite a colleague to see the loupes for themselves.

MAGNIFICATION

2.5× to 5.0×

WORKING DISTANCE

25–70 cm

PUPILLARY DISTANCE

52–78 mm

ILLUMINATION

Wired or wireless

WARRANTY

Lifetime

