




PRECISION LOUPES, FITTED TO YOU IN THE UK

See every detail Keep your posture


Deflection prism optics let you work looking straight ahead, on a frame adjusted to your eyes, your working distance and the way you sit or stand at work.

Look straight ahead, not down

Conventional loupes converge your eyes toward a single close point, and to see the field you tip your head down and hold it there. PENTAX Loupes are built around a different idea: your eyes stay level, and the optics do the bending.

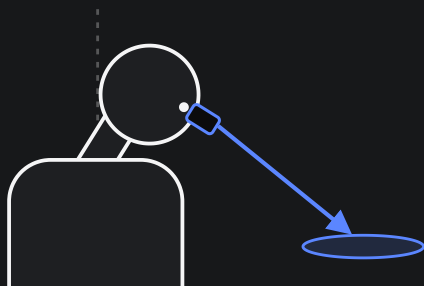

Conventional loupes. The head tips forward so the eyes can converge on one close point, and stays there for the whole procedure.

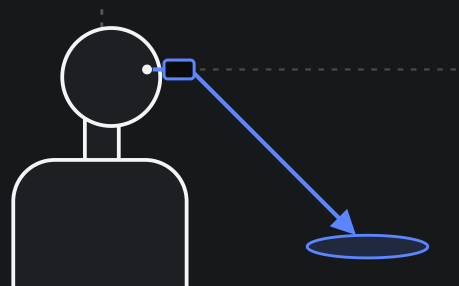

PENTAX Loupes. The head stays level and the eyes look straight ahead. Prisms inside the barrels bend the optical path down to the field.

Deflection prisms

Prisms inside each barrel redirect the optical path, much like working with a microscope. You look straight ahead at your work instead of angling your eyes down, which helps you hold a natural, upright posture through long procedures.

Parallel viewing

Conventional loupes converge your eyes toward a single close point. PENTAX Loupes are built for parallel viewing instead, designed to reduce eye convergence and help support an upright working posture.



"With PENTAX loupes, the posture of my neck really changed. I can look straight ahead, just like working with a microscope, and I enjoy parallel viewing with a very large field of view."

Dr Robotti · President, The Rhinoplasty Society of Europe

Built around you

Every measurement that matters is adjustable. It is set once at your fitting, then revisited as your eyesight or working posture changes.



25–70 cm

Working distance

Set to match your height, your chair and the way you sit or stand at work.



52–78 mm

Pupillary distance

Matched to your own eyes, not a factory average.



Height & angle

Line of sight

Position the barrels for a natural line of sight, then lock them in.



Optional lenses

Protective & prescription

Wear your glasses underneath, fit a prescription insert, or use the frame on its own.

Adjustable later, too. Pupillary distance and working distance are user-adjustable, so the fit can be revisited as your eyesight or working posture changes, without returning the loupes to the factory. Magnification is the one setting that is fixed: it is set when the loupes are built, which gives a consistent, predictable image.

Comfortable through long procedures

Posture is only part of it. The frame, the optics and the flip-up barrels are designed so the loupes stay comfortable and useful from the first case to the last.



Balance point

A lightweight frame with balanced weight distribution, so the loupes sit evenly and stay comfortable through long procedures rather than dragging on your nose.



A wide field of view

Japanese optics give a sharp, distortion-free image across a wide field, with the depth perception fine work needs. At 2.5× the field of view is 141 mm.



Flip up, flip down

Flip the barrels out of the way to switch instantly between magnified and normal vision, without taking the frame off.



"With the PENTAX Loupes I have a wide field of vision and they are light to wear. You can keep your focus while looking straight ahead, which minimises strain on my neck."

Dr Dhooghe · Surgeon



MAGNIFICATION

Choose your magnification

Five fixed magnifications, each with its own field of view. The choice comes down to what matters more for your work, a wider field or finer detail, and your specialist helps you decide at your fitting.

MAGNIFICATIONFIELD OF VIEW GUIDANCE

2.5×	141 mm	For those who favour field of view and depth of vision: general surgeons, urologists and hygienists.
3.0×	120 mm	For general dentists and ENT surgeons.
3.5×	106 mm	For anyone who wants a good middle ground between field of view and magnification.
4.0×	98 mm	For those who need good attention to detail while maintaining a decent field of view.
5.0×	80 mm	For those working on very small, fine details, such as endodontics or microsurgery.

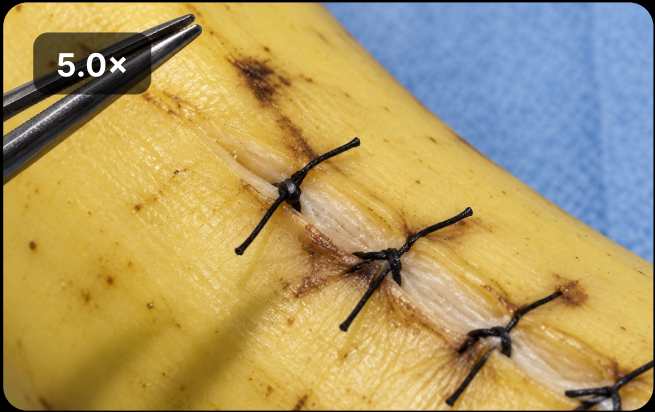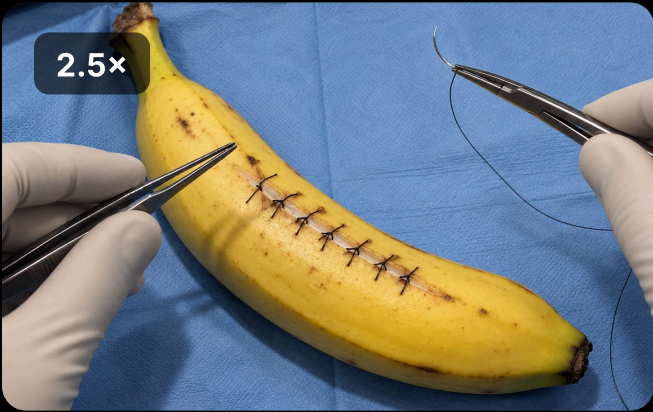

*These are estimates, results per person may vary.

2.5x–3.5x

Dentistry

Restorative and hygiene work at 2.5x to 3.5x. Aesthetic dentists often go to 5.0x.

3.5x–5.0x

Surgery

Fine, delicate procedures at 3.5x and above, across surgical specialities.

2.5x–3.5x

Veterinary

General and referral practice, depending on procedure and patient size.

Made for your eyes

Magnification is fixed and does not correct your everyday vision. There are three ways to pair PENTAX Loupes with your prescription.

Adjustable frame, no glasses

Many clinicians do not need separate glasses at all. The frame adjusts directly to your vision.

Your own glasses underneath

A different nose pad creates room to wear your existing glasses under the loupes.

Custom prescription insert

A clip-in insert carries your correction inside the PENTAX frame itself.



ILLUMINATION

Light exactly where you are looking

An optional LED headlight clips to the loupes and lights the field directly ahead of you.



Wired

Constant power from a backpack, for the longest lists.



Wireless

Nothing to clip to your clothing, for complete freedom of movement.

Already own a light you like? Brackets are available for other LED brands. Ask us about compatibility.

COLOUR

Four colours, pick yours



Silvergold



Blue



Black



Red

UK SERVICE

From enquiry to fitting

01

Tell us what you need

Reply to this email, or build your loupes online. It takes about a minute.

02

Talk it through

A specialist contacts you to discuss magnification, working distance and fit.

03

Demonstration & fitting

We arrange a demonstration and fitting with you in the UK.

AFTER YOUR FITTING

Every pair comes with a passport

A private page and a QR card in the case, holding your setup, your care guidance and a one-tap line to the UK team. Magnification, working distance and pupillary distance are recorded as they are set at your fitting, so any later adjustment starts from what was actually done.

Scan the card with any phone. No app, no account, no password. Use it for a quick check-in, a support request, or to invite a colleague to see the loupes for themselves.



WHO IT'S FOR

Dentists, prosthodontists and dental hygienists, surgeons across specialities, veterinary surgeons, and other professionals whose work requires enhanced close-up vision.

The best way to understand PENTAX Loupes is to wear them

Request a demonstration and fitting in the UK. A specialist will talk through magnification, working distance and fit with you, and prepare a personalised quotation for your configuration.

WEBSITE

pentaxloupes.co.uk

EMAIL

info@pentaxloupes.com

MAGNIFICATION

2.5× to 5.0×

WORKING DISTANCE

25–70 cm

PUPILLARY DISTANCE

52–78 mm

ILLUMINATION

Wired or wireless LED

